


**DERMASCOPE compact**  
contact/non-contact polarized dermoscopy







two color options available

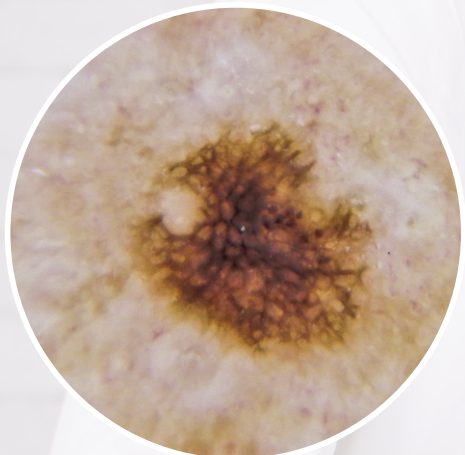

crystal clear images



safety wrist strap - standard accessory

## DERMASCOPE Compact handheld optical dermatoscope

DERMASCOPE Compact is a diagnostic device for optical, manual skin inspection. It utilizes contact or noncontact polarized dermoscopy method.

16 white LEDs as light source and three-way optical system comprised of glass magnifying lens, light polarizer and antireflex optics make this device ideal for patient screening. For greater mobility and autonomy rechargeable Li-Ion battery 600mAh is installed and can be charged over installed mini USB connector.

### Main advantages

- optical magnification cca 10X
- light intensity digitally adjustable in 3 levels
- non-contact and contact mode of operation
- polarizer for clearer and deeper images
- adapter for digital cameras
- all optics made of optical glass BK7 (no plastic lenses inside)
- color temperature 4500 - 5650K
- max light intensity 32500mCd
- CRI (Color Rendering Index) min 70
- compact dimensions 130x60mm
- two years limited warranty

### PROXIMA - medical technology

[www.proxima-medical.rs](http://www.proxima-medical.rs)

[office@proxima-medical.rs](mailto:office@proxima-medical.rs)

FB: ProximaMedTech IG: proxima\_medtech

phone/fax: +381 18 4291290, +381 18 4267469  
Moravska 7, Cokot 18250 Novo Selo, Nis, Serbia



16 white LEDs for best visibility



digital dermoscopy system



adapter for digital camera