

RWSU 70

4MHz Radio Wave Surgical Unit

powered by

proxima
medical technology



SA SAMO TRI MOĆNA MODA RADA ZA SVE VRSTE INTERVENCIJA:

SOFT CUT - čisto rezanje sa minimalnim oštećenjem tkiva

STRONG CUT - snažno rezanje sa izraženijim efektima koagulacije

SOFT COAG - soft kontaktna koagulacijat

WITH ONLY THREE POWERFUL MODES FOR ALL INTERVENTIONS:

SOFT CUT - clean cutting with minimal tissue damage

STRONG CUT - powerful cutting with stronger coagulation effect

SOFT COAG - soft contact coagulation



GENIJALNA JEDNOSTAVNOST KORISNIČKOG

- Selektor za izbor moda rada sa svetlosnom indikacijom u boji (žuto ili plavo)
- Povećanje ili smanjenje snage samo okretom dugmeta sa pozadinskim svetlom u boji
- Aktivacija snage samo nožnom komandom za veću manuelnu preciznost u radu

INGENIOUS SIMPLICITY OF USER INTERFACE:

- Mode selector knob for operating mode with color coded light indication (yellow or blue)
- Increase or decrease output power with just a turn of a button with color backlight
- Power activation by footswitch only for higher manual precision

SIGURNOST NA PRVOM MESTU:

- Implementirano sigurnosno kolo neutralne elektrode sa audio i svetlosnom signalizacijom i prekidom izlazne snage do uspostave regularnog kola neutralne elektrode

Korišćenje samolepljivih neutralnih elektroda


- smanjuje mogućnost kros-kontaminacije kod pacijenata

SAFETY FIRST:

- Implemented neutral electrode safety circuit with auditive and visual signalization and output power interruption until regular neutral electrode safety circuit is established

Usage of self-adhesive neutral electrodes

- lowers the possibility of cross-contamination between patients

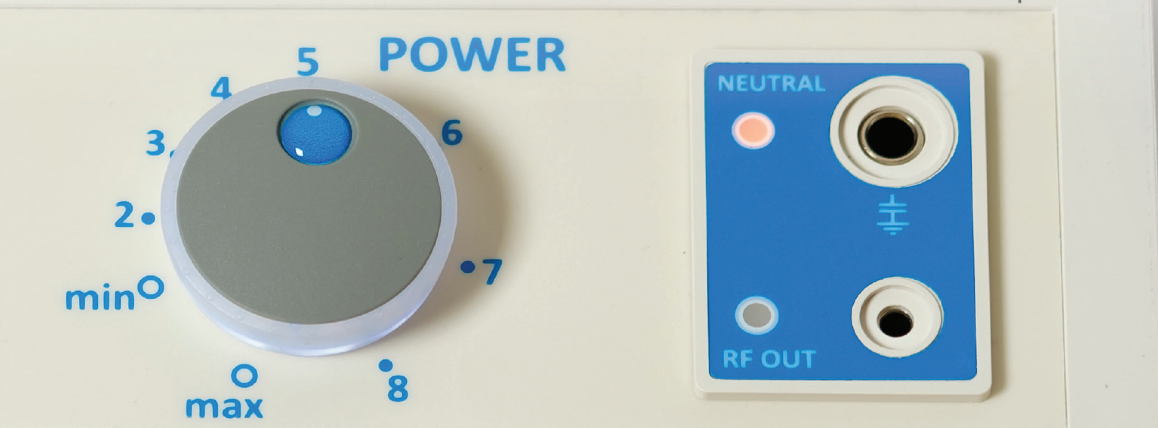

PRIMENJENO INŽENJERSTVO ZA LAKO KORIŠĆENJE:

- Konektori za aktivnu i neutralnu elektrodu na prednjoj strani, konektori za nožni prekidač i napajanje na zadnjoj strani - omogućuju čistiju radnu površinu bez umotavanja kablova

Male dimenzije i desktop dizajn čine RWSU 70 pogodnim kako za radne površine tako i za odgovarajuća kolica i lako korišćenje između više specijalnosti (npr. dermatologija i ginekologija u istoj klinici)

- Male dimenzije i desktop dizajn čine RWSU 70 pogodnim kako za radne površine tako i za odgovarajuća kolica i lako korišćenje između više specijalnosti (npr. dermatologija i ginekologija u istoj klinici)

APPLIED ENGINEERING FOR USABILITY:

- Connections for active and neutral electrode on the front, connections for footswitch and power supply on the back - enabling clean working space and no messy cables

Small dimensions and desktop design make RWSU 70 suitable for desktops as well as

- appropriate trolleys and easy usage between several specialties (eg. dermatology and gynecology at the same clinic)

PRIBOR

Za uređaj RWSU 70 dostupan je širok spektar monopolarnih elektroda za višekratnu upotrebu:

Kratke (cca 7cm) monopolarne elektrode (za dermatologiju, ginekologiju, plastičnu hirurgiju, maksilofacijalnu hirurgiju, oralnu hirurgiju, ORL...)

Dugačke (cca 12cm) monopolarne elektrode (za ginekologiju, plastičnu hirurgiju, ORL...)

Sve elektrode su za višekratnu upotrebu i mogu se sterilisati u autoklavu (parnom sterilizatoru).


ACCESSORIES

For RWSU 70 device a wide range of monopolar electrodes for multiple use is available:


Short (cca 7cm) monopolar electrodes (for dermatology, gynecology, plastic surgery, maxillofacial surgery, oral surgery, ENT...)

Long (cca 12cm) monopolar electrodes (for gynecology, plastic surgery, ENT...)

All electrodes are for multiple use and can be sterilized in autoclaves (hot steam sterilizers).

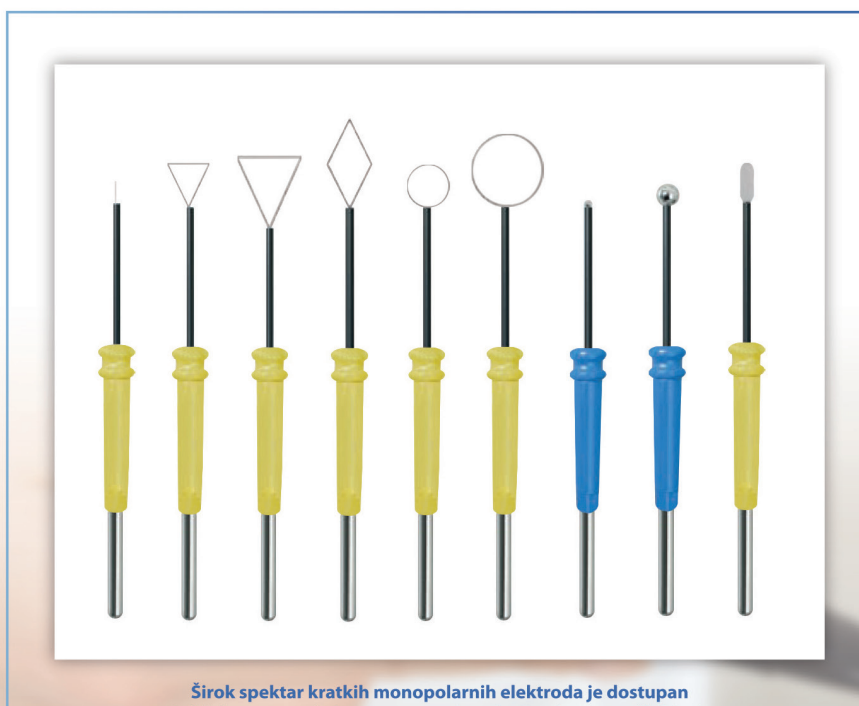

Širok spektar kratkih monopolarnih elektroda je dostupan
Wide range of short monopolar electrodes is available

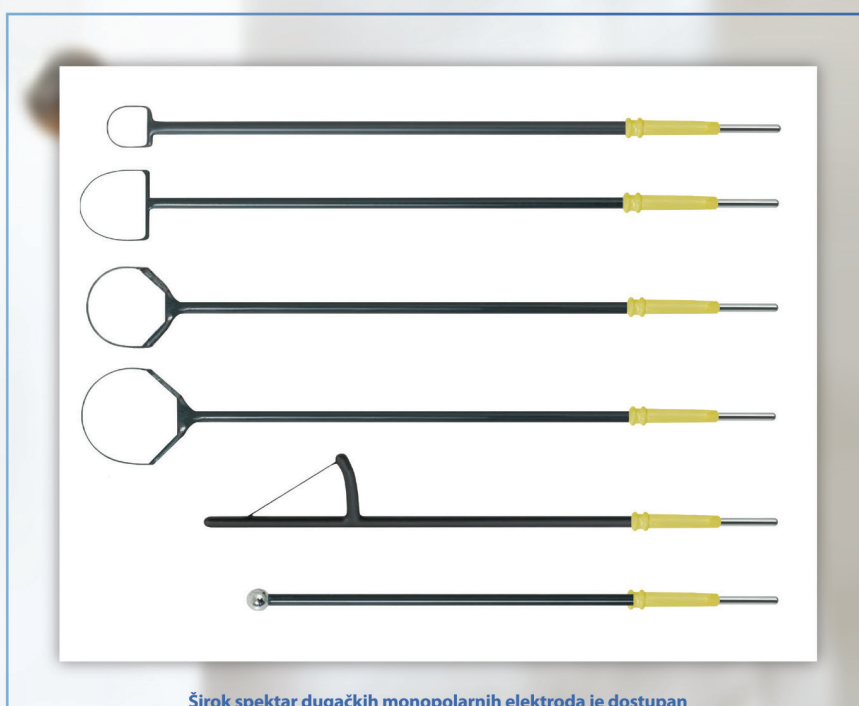

Širok spektar dugačkih monopolarnih elektroda je dostupan
Wide range of long monopolar electrodes is available

TEHNIČKE KARAKTERISTIKE

TECHNICAL CHARACTERISTICS

NAPAJANJE POWER SUPPLY	230V 50Hz
MODOVI RADA MODES OF OPERATION	Soft Cut Strong Cut Soft Coag
IZLAZNA SNAGA OUTPUT POWER	70W RMS @ 4MHz, sinusni oblik / sine wave
RADNA FREKVENCIJA OPERATING FREQUENCY	4MHz
SIGURONOSNO KOLO NEUTRALNE ELEKTRODE SAFETY CIRCUIT OF NEUTRAL ELECTRODE	Da / Yes

ŠIFRE ZA PORUČIVANJE

ORDER CODES


UREĐAJ RWSU 70 SA KOMPLETNIM STANDARDNIM SETOM AKSESORIJA RWSU 70 DEVICE WITH COMPLETE STANDARD SET OF ACCESSORIES	PX-AP-2051
KOLICA SA TOČKOVIMA ZA RWSU 70 TROLLEY WITH WHEELS FOR RWSU 70	PX-AP-2051-C1





PROXIMA - medical technology

 office@proxima-medical.rs

 +381 18 4 291 290, +381 18 4 291 288
Fax: +381 18 4 267 469

 Moravska 7, Čokot
18250 Novo Selo, Niš, Republic of Serbia



1304



 [instagram.com/proxima_meditech](https://www.instagram.com/proxima_meditech)

 [facebook.com/ProximaMedTech](https://www.facebook.com/ProximaMedTech)

Distributor: