

TEHNIČKE KARAKTERISTIKE ES200RKB

Napajanje	230V 50Hz
Dimenzije	309x290x129mm (DxŠxV)
Max srednja snaga	monopolarni rez 200W na 500Ω monopolarna koag. 150W na 500Ω bipolarna koag 120W na 200Ω
Izlazna frekvencija	490kHz sinusni signal

TECHNICAL DATA ES200RKB

Power supply	230V 50Hz
Dimensions	309x290x129mm LxWxH
Max average power	monopolar cut 200W on 500Ω monopolar coag 150W on 500Ω bipolar coag 120W on 200Ω load
Output frequency	490kHz sine wave



Ergonomski kolica za aparat i pribor
Ergonomic trolley for device and for accessories

ŠIFRE ZA PORUČIVANJE

ES200RKB aparat	PX-AP-2028
Neutralna silikonska elektroda za odrasle	PX-AP-1001-06
Neutralna silikonska elektroda za decu	PX-AP-1007-06A
Kabl za neutralnu silikonsku elektrodu	PX-AP-1001-07
Nožni prekidač	PX-AP-1001-01
Nožni prekidač - economy	PX-AP-1001-01A
Mrežni kabl	PX-AP-11
Kolica za aparete i pribor	PX-AP-2027-01

ORDER CODES

ES200RKB device	PX-AP-2028
Neutral silicon electrode for adults	PX-AP-1001-06
Neutral silicon electrode for children	PX-AP-1007-06A
Cable for neutral silicon electrode	PX-AP-1001-07
Double footswitch	PX-AP-1001-01
Double footswitch - economy	PX-AP-1001-01A
Mains cable	PX-AP-11
Trolley for device and accessories	PX-AP-2027-01

PROXIMA d.o.o.
Moravska 7, Čokot 18250 Niš
Tel: ++381 18/ 4291290 ; 4291288
Tel/fax: ++381 18/ 4267469

www.proxima-medical.rs
office@proxima-medical.rs



elektrohirurški generator
ES200RKB

HF ELECTROSURGERY

KARATKTERISTIKE ES200RKB / ES200RKB FEATURES

- rez, koagulacija i bipolarna koagulacija
- aktiviranje ručnim operacionim držaćem ili nožnim prekidačem
- siguronosno kolo neutralne elektrode i alarmiranje sa prekidom svih funkcija
- max izlazna snaga 200W u monopolarnom i 120W u bipolarnom režimu rada
- jednostavno i precizno zadavanje snage obrtnim potenciometrima i izbor vrste rada
- monopolarno-bipolarno putem tastera na folijskoj tastaturi
- namenjen za ambulantne i uslove kliničkih operacionih sala
- širok spektar pribora omogućava veoma veliki broj aplikacija
- efikasno rezanje, koagulacija i bipolarna koagulacija i u otvorenoj i u endoskopskoj hirurgiji
- kompaktan dizajn i povezivanje svih vrsta pribora odjednom na prednjoj strani aparata, sa zasebnim priključcima za monopolarni i bipolarni

- cut, coagulation and bipolar coagulation
- hand operational switch or footswitch controlled
- neutral electrode safety circuit and alarm; function break in case of alarm
- max output average power of 200W in monopolar and 120W in bipolar operation
- simple and precise power adjustment by rotary potentiometers and mode choosing between monopolar and bipolar by taster on the front panel
- application in ambulatory conditions and ORs
- wide range of HF accessories makes possible great variety of applications
- efficient cutting, coagulation and bipolar coagulation both in open and endoscopic surgeries
- compact design and connection of all accessories on the front panel of device, separate plugins for monopolar and bipolar

dvostruki nožni prekidač - „high value“
 double footswitch - „high value“



dvostruki nožni prekidač - „high value“, vodoootporan,
 sa prekidačima za Rez i Koagulaciju prema EN60601-2-2
 double footswitch - „high value“, waterproof,
 with switches for Cut and Coagulation acc EN60601-2-2



* izrađen od aluminijuma i plastike

neutralna elektroda za odrasle sa kablom
 adult neutral electrode with cable



HEBU medical, Tutlingen - Nemačka

PROXIMA d.o.o. je posvećena kvalitetu svojih proizvoda, kako elektrohirurških jedinica, tako i elektrohirurškog pribora kao osnovnog alata koji dolazi u interakciju i sa pacijentom i sa hirurgom. Upravo zahvaljujući našem nastojanju da klijentima pružimo kvalitet postali smo zvanični distributeri jednog od najvećih svetskih proizvođača elektrohirurškog pribora - HEBU medical. Sav pribor je izrađen vrhunskom tehnologijom u skladu sa evropskim i svetskim standardima i opremljen CE i FDA sertifikatima. U ponudi je čitav spekter pribora za otvorenu i za laparoskopsku i endoskopsku elektrohirurgiju.




Tropinski Valleylab konektor



Proxima/Martin/Berchtold
konektor



Erbe ACC/ICC konektor

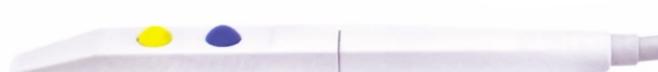

Ručni operacioni držači za 400 sterilizacija



Veliki izbor aktivnih elektroda
 (više od 100 različitih elektroda u ponudi)

HEBU
Non-Stick



Širok spekter bipolarnih pinceta u ponudi



* PROXIMA je distributer pribora za elektrohirurgiju renomiranog nemačkog proizvođača HEBU medical